

permobil



NEW KID IN TOWN

The all-wheel suspension system of the K450 MX provides a smooth ride, whether cruising the halls at school or heading out to the playground. The power seat-to-floor function lets children get down on the ground, perfect for story time, transfers and playing with others. The seat-to-floor height ranges from 3" to 26," making it easy to sit at a table or counter and even to reach things down low. Plus, the K450 MX features R-Net electronics and has an array of power functions, including: seat elevator, seat-to-floor and optional power tilt. The modular seating system has center mount adjustable legrests and comes with a free growth kit—something every kid needs.


K450
MX



POWER FUNCTIONS



SHROUD COLORS

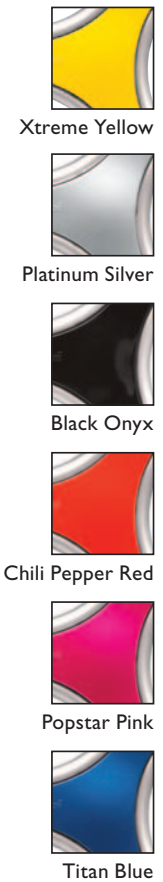

R-Net electronics for improved programmability



Optional 45° power tilt from any seat position



Tool-free adjustable legrest




All-wheel suspension system for increased balance, improved ride and safety



Large flat free tires for improved performance and reduced maintenance



Uni-track system allows for extensive customization

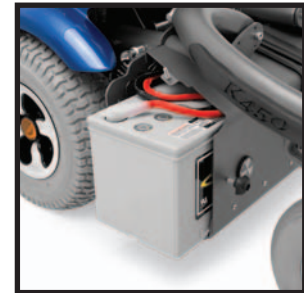

Easy battery access

SPECIFICATIONS

Maximum Speed	5.2 MPH	Power Seat Elevator	8"
Max User Weight	125 lbs	Seat Widths	10", 12", 14", 16"
Ground Clearance	3"	Seat Depths	10" - 12", 12" - 14", 14" - 16", 16" - 18"
Battery Capacity	Group 34	Backrest Heights	14" - 20"
Minimum Width	24.5"	WC-19 Crash Test Approved	Yes
Base Length	33"	Seat To Floor Height w/ Tilt	14" - 18"
Weight of Wheelchair w/Batteries	311 lbs	Tilt Angle	45° Power
Uphill Angle	18°	Electronics	R-Net
Downhill Angle	18°	Alternative Drive Controls	Yes w/R-Net
Side Slope	16°	Effect of Driving When Tilted	Low speed at 20°



THE POWER OF MOBILITY

Permobil Inc.
300 Duke Drive
Lebanon, Tennessee 37090
800.736.0925
www.permobil.com