



## DRIVE IN STYLE

The Permolock C3 allows you to drive or ride while sitting in your Permobil. This device secures your wheelchair in place inside the vehicle and is designed for use with the Permobil C300 and C350 wheelchair bases. You can adjust both the angle and height of the Permolock C3, which helps when mounting it to a vehicle with a sloped driver's position. A locking device mounted to the wheelchair chassis and a locking plate that is mounted to the vehicle ensure the wheelchair does not break free during a collision or because of some other sudden movement. Plus, you can either pull forward or back into the locking plate.

**PERMOLOCK**  
C3



### SPRING-LOADED PINS

Two fixed pins attach to the Permlock C3 base to allow for a safe and secure ride. When an obstacle is encountered, they push up into the chassis.



### LOCKING PLATE

Easy-to-install docking station that holds the chair firmly and securely in position.



### ADJUSTABLE BASE

The height of the base can be adjusted in 5mm increments, and the angle adjustment utilizes a spring to help with mounting it on a slope.




### VEHICLE CONTROL BUTTON

Located on the dash, it uses three lights—red, amber and green—to indicate if the chair is locked into place. Simply press and hold to unlock the wheelchair.



### EMERGENCY RELEASE


Locking mechanism can be opened manually with the emergency handle on the front edge of the locking plate.

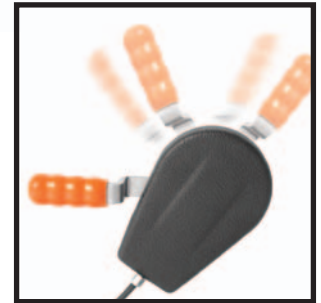

### MANUAL RELEASE

Wheelchair can be released manually in the event of vehicle power failure.

## ACCESSORIES

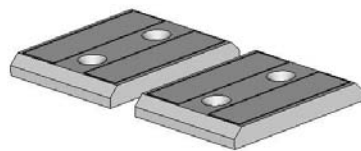

### INSTALLATION SPACER

Used for installing Permlock on a sloping surface.



THE POWER OF MOBILITY

Permobil Inc.  
300 Duke Drive  
Lebanon, Tennessee 37090  
800.736.0925  
www.permobil.com

\* Please consult owner's manual before using Permlock.

Permlock must be installed by a certified NMEDA representative.